Tecnorita
by **BH**

Instrucciones de uso y garantía
Operating instructions and guarantee
Mode d'emploi et garantie
Instruções de utilização e garantia
Bedienungsanleitung und Garantie
Manuale d'uso e garanzia

BIENVENIDA

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WELCOME

We would like to thank you for purchasing our product and for the trust you have placed in Tecnovita. All of our products have been tested and submitted to strict quality controls. Please read this manual before using the product.

BENEFITS

- THE BEST TECHNOLOGY TO ELIMINATE CELLULITIS AND ORANGE PEEL SKIN;
- HIGHLY EFFECTIVE MUSCLES TONING;
- IMPROVES BLOOD CIRCULATION AND ACCELARATES METHABOLISM;
- SCULPTS THE BODY;
- IMPROVES THE ABSORTION OF SKIN PRODUCTS.

ULTRASONIC USE**How ultrasonic works?**

The high frequency ultrasonic wave technology helps to improve the cell membrane metabolism, contributing to better blood circulation This technology increases the oxygenation of skin cells providing a completely safe and painless way to combat the sagging skin, to remove cellulitis and to tone muscles. The ionic technology improves facial cleansing, helping to remove grease, dirt and traces of makeup.

What are the advantages over other methods?

1MHz ultrasonic waves are safe, painless and non-invasive and they aid sonophoresis, allowing the nutrients of skin care creams to penetrate deeper.

Is it painful?

It is painless and safe.

ELECTROSTIMULATION METHOD

EMS, electric muscular stimulation is a widely recognized method that it has been used on the sports field for a long time. It allows enjoying the exercise and increasing muscle capacity by developing strength and improving tone effortlessly. Thanks to the DualCare design it can tone muscles anytime, anywhere.

TECHNICAL CHARACTERISTICS

- 3- in-one body treatment;
- Body ultrasonic with 1MHz frequency;
- Body EMS system for body toning of muscles and sculpts the body;
- Infrareds that help improving blood circulation and metabolism;
- Includes power adapter.

SPECIFICATIONS

Output power: 1 MHZ

Weight: 140 gr.

Power adapter. 100-240 VAC 50-60 Hz 0.2A max

Dimensions : L x W x H (190 x 80 x 60 mm).

RECOMMENDATIONS FOR USE

Precautions: Consult your doctor before start using this appliance. This warning is especially important for people with health problems. Use by those under 18 years of age is not recommended. Keep this appliance out of children's reach.

- This appliance is designed for home use;
- This appliance is design to be used in the ultrasonic mode in arms, abdomen, flanks, gluteus, hips, thigh and calf muscle.

WHEN IS NOT RECOMMENDED THE USE OF DUAL CARE?

- if you suffer from hyperpigmentation, broken capillaries;
- if you have sunburn;
- if you suffer from severe skin problems;
- active skin infections;
- facial paralysis;
- If you are pregnant or think you might be;
- If you are breast-feeding;
- If you are within your menstruation period;
- Undergoing fertility treatments;
- Couperosis;
- Osteoporosis;
- Fibromyalgia;
- malignant tumors;
- diabetes;
- kidney problems;
- physical or mental deficiencies;
- metal implants or prostheses in the area of application;
- scarring in the area of application;
- Physical malformations;
- Epilepsy;
- long-term treatment with medications: anti-inflammatories;
- Anticoagulants, antibiotics or antihistamines;
- treatment for psychological or psychiatric disorders;
- heart diseases;
- metal implants;

We recommend to use the device on a clean skin, with no make-up in order to facilitate ultrasonic penetration.

In the ultrasonic mode always spread sufficient gel to avoid poor transmission or damage of the device.

To avoid skin burn, never stop moving the device onto the skin when functioning.

Do not use the device for more than 15 minutes in the ultra ultrasonic mode.

Do not use the device for more than 20 minutes in the LED mode.

Always use the safety glasses in the LED mode.

Do not use the device if the cable is damaged.

To unplug the device do not pull the cable.

Do not touch the device when your hands are wet.

Never submerge the device into water.

Always clean the device with a smooth cloth after each use.

Do not expose the device to high temperature or intense magnetism.


If you notice an incorrect functioning of the device, do not try to open, disassemble or modify it. No authorized repairs can damage the device and void warranty. Please contact TS.

PACKAGE CONTENT

1. DEVICE
2. ULTRASOUND GEL
3. POWER CABLE
4. ELECTRODES CABLE
5. ELECTRODES

DESCRIPTION

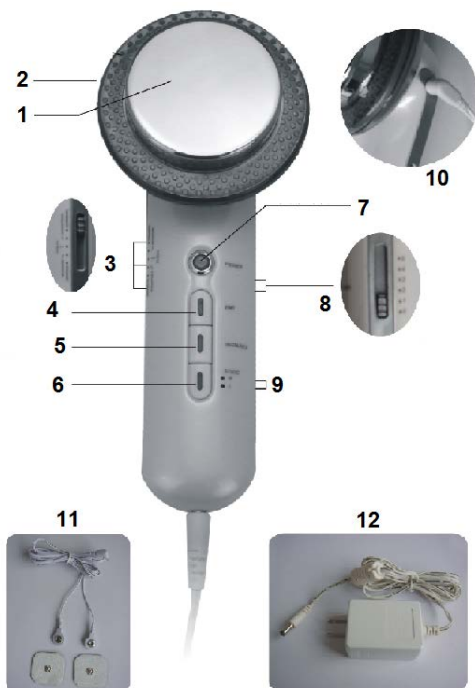1. Ultrasonic wave massage probe;
2. Far Infrared Light;
3. EMS mode selection;
4. EMS mode;
5. Infrared mode;
6. Ultrasonic mode;
7. Power button
8. EMS intensity
9. Ultrasonic intensity
10. EMS socket
11. Cable and contact pads
12. Power cable

QUICK USER GUIDE

EMS + Infrared + Ultrasonic >>>> Multiplied effect

1. Please Clean your skin before application;
2. Use EMS for 15 minutes on the are that needs slimming;
3. Apply abundantly the ultrasonic gel;
4. Use the device on ultrasonic and infrared mode until the gel is totally absorbed into the skin. (Max 15 minutes).

INSTRUCTION FOR USING ULTRASONIC

Step 1

Cleanse your skin before application

Step 2

Spread the body lotion evenly (apply the product separately for each section in order to avoid air-drying).

Step 3

Plug in the Dual Care device.

Step 4

Turn power ON and a red light will blink

Step 5

Press the SONIC button – once for low intensity, twice for high intensity. Use the Low intensity on the face; high intensity can be used for the rest of the body.

Step 6




Spread the gel on the area you wish to apply massage. Keep the skin moisturized.

Step 7

Follow instructions of images for each part of the body.


ARMS

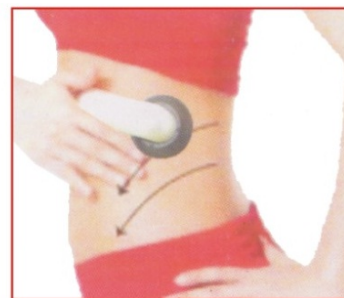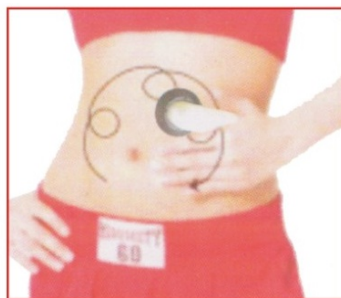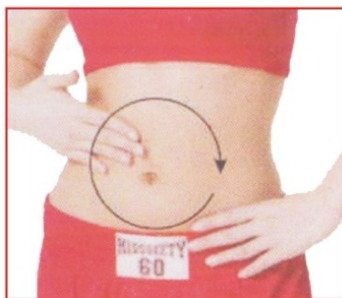WAIST AND ABDOMEN

HIPS AND GLUTEUS

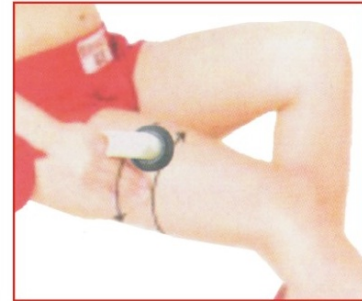LEGS AND FEET

For a correct and safe use of DualCare please consider the following application times:

Daily session: 15 minutes maximum

A session is intended as the use of DualCare in the same day on the several parts of the body.

A minimum of 48 hrs must be delayed between sessions: (i.e. Monday, Wednesday and Friday).

Every 24 sessions (one month and a half) it is recommend to give body a rest of another month and a half.

VERY IMPORTANT

*****Never exceed times of application (within 15 minutes in total daily).** Exceeding application times can cause fatigue and/or skin irritation.

*****Please observe the resting time of 48 hours between one session and another.** This is the time your skin needs to rest.

INSTRUCTIONS FOR USING EMS

1. Clean the skin area that will be treated;
2. Connect the EMS cable with electrodes. Place the electrodes on the area that will be treated. Positive and negative electrodes works simultaneously;
3. Select 0 intensity level;
4. Connect the EMS cable with the device
5. Press the ON/OFF button to switch on the device. A red light will lite;
6. Select the EMS mode (Tapping, Massage, Knead, Scrapping, Slimming)
7. Press the EMS button
8. Select the intensity of the massage (1-5)
9. The device will automatically switch off after 10 min. If you wish to continue, please turn on device and repeat the procedure.
10. After the treatment, adjust the mode selection key to the no-mode position and the 0-intensity position.

INSTRUCTIONS FOR USING INFRARED

When using the Far Infrared, it is important to apply skin care product on the skin. The skin product is both conduction medium and skin nutrition.

1. Plug the device to electric power;
2. Press the ON/OFF button to switch on the device. A red light will lite;
3. Select the Infrared mode;
4. Point the infrared on the areas that will be treated.
5. Turn off the device making sure that all functions are turned off.

FRECUENTLY ASKED QUESTIONS

How ultrasonic works?

The high frequency ultrasonic wave technology helps to improve the cell membrane metabolism, contributing to better blood circulation This technology increases the oxygenation of skin cells providing a completely safe and painless way to combat the sagging skin, to remove cellulitis and to tone muscles.

The ionic technology improves facial cleansing, helping to remove grease, dirt and traces of makeup.

What are the advantages over other methods?

1MHz ultrasonic waves are safe, painless and non-invasive and they aid sonophoresis, allowing the nutrients of skin care creams to penetrate deeper.

Is it painful?

It is painless and safe.

Is it normal to hear a slight buzzing in your ears while using DualCare?

Yes, some users may detect a slight buzzing sound during its application.

Can it be used on all parts of the body?

The 1 Mhz DualCare technology has been specifically designed for its use on arms, waist, hips and legs.

After how many sessions will I begin to see results?

You should notice results after the first five sessions, however results may vary depending on the skin of each individual.

How long should I wait between sessions?

It is very important to respect the 48 hour break between sessions, as this is the time that your skin needs to rest.

Are there any side effects?

If you are healthy, DualCare has no reported side effects, nevertheless in some people it may produce a reddening of the skin that can last for a few minutes.

What is ultrasound gel for? When is it used?

The ultrasound gel is essential, ALWAYS use it before beginning a session. This will ensure that the ultrasonic waves penetrate correctly and prevent damage to the unit.

Can I use DualCare without applying gel?

Never, because the gel acts as a conductor for the ultrasonic waves on your face. If you do not use gel then the DualCare application loses its effectiveness and the machine can be severely damaged.

Can I use any type of gel?

We recommend using only TECNOVITA gel, as the use of any other type may damage the machine.

Is it necessary to drink water?

No, it is not strictly necessary but it is advisable to hydrate yourself. The recommended quantity is at least 2 litres during treatment, as this helps to hydrate the skin.

Can I use DualCare if I am pregnant?

NEVER use DualCare if you are pregnant or suspect that you might be.

Can DualCare be used during lactation?

Do not use it if breastfeeding.

Can DualCare be used on the breasts?

Never apply it to the breasts.

Can DualCare be used with skin conditions?

No, It should not be used if you have any serious skin condition. Consult your doctor.

Can I use DualCare if I have hearing implants?

No, it should not be used if you have hearing implants fitted.

Can I use DualCare if I have implants in my mouth, gold teeth, metal implants or similar?

Do not use it if you have these types of implant, as it may cause discomfort or tenderness in the area of application.

Can I use DualCare when I am bathing or showering?

Under no circumstances should the handset or the machine come into contact with water due to the risk of electrocution.

How do I clean DualCare?

Never immerse DualCare in water to clean it. Use a cloth slightly dampened with water.

Warranty

We suggest you to read this document carefully in order to fully understand the cover that the manufacturer offers as well as your responsibilities as owner.

GUARANTEES

GENERAL INFORMATION

EXERCYCLE, S.L. would like to thank you for purchasing this product and the trust you have placed in us. Our products have been subjected to rigorous quality checks. Nevertheless, should you experience any problem with this appliance during the guarantee period then please read these warranty conditions and then get in touch with your local authorized dealer.

This guarantee certificate has been drafted in accordance with that stated in Directive 1999/44/EC of the European Parliament and of the Council 25th May, 1999, regarding certain aspects of the sale of consumer goods and associated guarantees, incorporated into Spanish Law by virtue of Law 23/2003, 10th July; Guarantees on the sale of consumer goods.

To this effect, EXERCYCLE, S.L. VAT nº B-01029032, with business address c/ Zurrupitieta 22, Vitoria, Spain, guarantees this product under the following terms and conditions:

1. GUARANTEE PERIOD

The product covered by this guarantee certificate is protected for a period of TWO (2) YEARS as from the date of receiving the product, which will be taken as the date that appears on the sales receipt for the product or on the corresponding delivery docket should this be later.

It is essential that this be presented in the event of possible repairs. The policy regarding the warranty period may vary between different countries.

The stated warranty period is only applicable to Europe.

2. RIGHTS UNDER THE GUARANTEE

- Repair or replacement of the product which, in such instance, will be carried out within a reasonable period and without significant inconvenience to the consumer, taking into account of the nature of the goods and the purpose for which the consumer required them.
- Those cases in which the repair or replacement represents a disproportionate cost for the seller in comparison with other forms of remedy will be exempt assuming these do not prejudice the consumer.
- Repairs will be free for the consumer with regard to labour, materials and carriage.
- In the event that the machine is repaired or replaced with a new one, the guarantee period will be suspended until the new appliance has been delivered and will carry at least six months of guarantee, however under no circumstances will a new 2-year period be initiated.

The buyer may, as is apparent, require (i) a reduction in price or proceed to a (ii) recession of the contract when the repair or replacement is not carried out within a reasonable time or with significant inconvenience to the buyer. The reduction in price (i) will be proportional to the difference between the value that the goods would have had if they had been delivered in accordance with the contract and the actual value that the goods were worth at the time of the said delivery. Recession of the contract (ii) will not proceed when the non-conformity is of minor importance.

3. THE GUARANTEE DOES NOT COVER

- Defects attributable to use other than that established in the User Manual.
- Professional use or commercial leasing or hiring of the product.
- Poor installation, poor assembly or breakage of components due to improper installation, (the User Manual gives detailed instructions on how to install and use the machine).

- Improper use of the product (use without applying gel, exceeding operating times, exposure to high moisture levels, etc.) or problems arising from lack of maintenance (accumulation of dirt on any part of the product or its accessories, etc.). (The User Manual provides advice on how to install and use the product).
- Repairs carried out by persons not authorized by EXERCYCLE, S.L. or appliances that have been opened or tampered with by persons not authorized by the manufacturer.
- Deviations from the standard product, such as light vibrations, noises or breakage of components that do not affect the performance of the machine and which cannot be attributed to manufacturing faults.
- Scratches, scrapes, abrasions, knocks falls, etc. caused by carelessness.
- Defects on the appliance caused by damage during transit.

4. CLAIMING AGAINST THE GUARANTEE

- Claims made for nonconformities will be limited to three years as from the date that the goods were delivered.
- The consumer must report the existence of the defect within a period of 2 months as from the date that the defect was discovered.
- This guarantee does not affect the consumer's statutory rights under any applicable law, unless it is incompatible with the application of other actions derived from compensation sought for hidden faults in the contract of sale.
- With regard to queries that may arise from the interpretation or application of this guarantee, these will be governed by the corresponding national law.

Any discrepancy arising from the interpretation, validity or application of this contract will be subject to the courts and tribunals of Vitoria, Spain.

Certificates

By virtue of this document, Exercycle S.L. declares that this product meets the essential requirements and other relevant provisions of Directives 2004/108/EC and 2006/95/EC

Please keep the receipt as proof of purchase.

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BH China Co., Ltd.